

Elliptic Reusable Bimanual I/A System



Product Features

- Titanium lightweight, colour-coded Handpieces
- Standard male/female Luer connector
- Compatible with Phaco or I/A System
- Smooth bullet shaped tip
- Rough tips available for capsule polishing
- Swiss Design / Swiss Made

Based on 21G, elliptic cannula system

- major axis outer diameter 0.94 mm, inner diameter 0.74 mm
- minor axis outer diameter 0.62 mm, inner diameter 0.42 mm

Irrigation

BIO 128/T

1 Irrigation Handpiece with curved cannula
Proximal TWIN-PORT® with diameter of 0.55 mm
per port



Aspiration

BAO 110/T

Aspiration Handpiece with curved cannula,
elliptic port 0.30 mm, roughened tip



Details



Important User information

The BIO 128/T and the BAO 110/T-system of Bürki are an elliptic I-/A-system. The higher impermeability due to the elliptic canula forms in the paracentesis leads to substantial less efflux. In consequence the inner eye pressure can increase substantially.

To prevent any damage to the eye, a reduced height of the irrigation liquid bottle is strongly recommended and the pressure conditions must stay in the permanent focus of observation and, if necessary, adjustment. Nevertheless it's important to prevent any collapse of the eye.

Version 06/2015

About Bürki
Bürki plans, develops, produces and certifies
microsurgical instruments in SWISS QUALITY
DIN EN ISO 13485

Bürki inno med AG
Industriestrasse 67
9443 Widnau
Switzerland

Tel +41 (0)71 744 88 66
Fax +41 (0)71 744 88 67
info@buerki.eu
www.buerki.eu

CONTATTI DISTRIBUTORE

SIR Oftalmica Srl

Via Torriani 19/C

22100 - COMO

Tel. : 031-570869

info@siroftalmica.com


www.siroftalmica.com

P.IVA : 01814920136

Contatti specialists:

A. Molteni: 328.1808176

L. Ronchi: 335.6167324

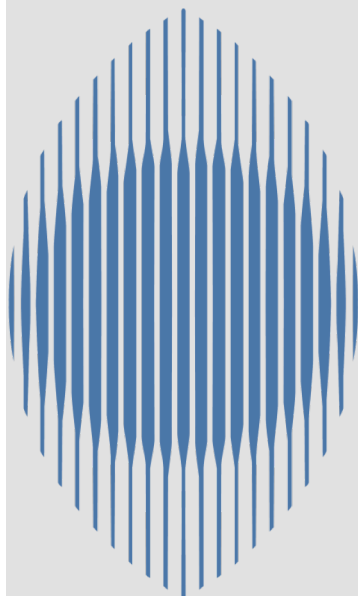

SIR
oftalmica

FERRI CHIRURGICI – PRODOTTI MONOUSO
APPARECCHIATURE DI DIAGNOSTICA E TRATTAMENTO
LETTINI E SEDIE CHIRURGICHE - RIUNITI